



# ADSEAL Plus

## Root Canal Sealer



Experience  
the Upgraded Performance!

ADSEAL Plus is more economical and convenient  
Resin based-sealer, which encompasses all advantages of ADSEAL.



### Economical Improvement with 2 Types of Tips


ADSEAL Plus offers a **green-colored dispense tip** that enhances control of achieving the desired amount of sealer and an auto-mix tip (optional, sold separately), which provides users with easy and convenient mixing.



### Excellent Physical Properties


#### Exceptional Flowability

ADSEAL Plus provides **the most appropriate flowability** for RCT when base and catalyst are mixed.

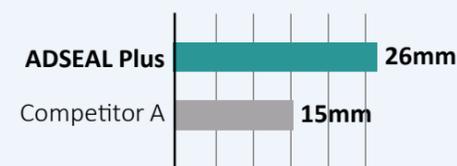


#### Accurate Diagnosis

ADSEAL Plus shows about **25% higher radiopaque** than competitor A.

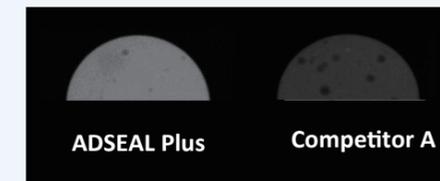

### Stability & User Friendly Grip

ADSEAL Plus syringe has a finger-friendly design with a reduced strain of 24mm flange size, while maintaining a **same weight of ADSEAL (13.5g)**

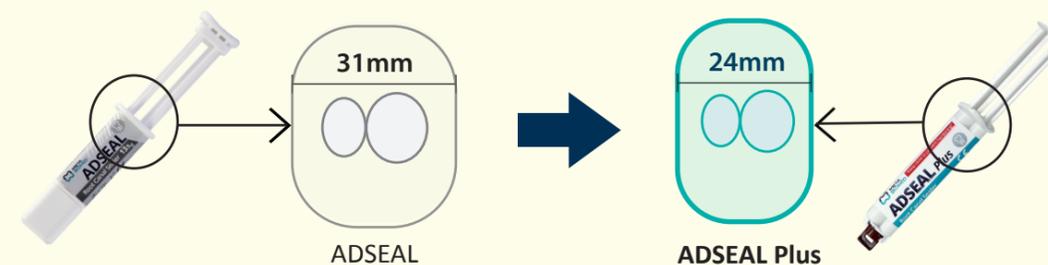

### Unparalleled Aesthetic Enhancement

ADSEAL Plus contains **Zirconium Oxide** which doesn't cause tooth discoloration. (**Bismuth-free**)

### High Biocompatibility & Increased Safety

ADSEAL Plus contains Calcium Phosphate, it offers **biocompatibility and safety** for humans.

#### SPECIFICATION

- Net weight : 13.5g(Base 9g, Catalyst 4.5g)
- Setting time : 50 min
- Radiopacity : 4.4mm(Al)
- Exp Date : 2years
- Auto mix & endo tips (sold separately) : 20EA/1pkg

#### ITEM CODES

- ADSEAL Plus: 303002
- ADSEAL Plus pre mixed tip: 303002TIP

#### COMPONENTS

- 13.5g Dual syringe(Base 9g, Catalyst 4.5g)
- Mixing pad
- Spatula
- Dispense Tip& Cap
- ADSEAL Plus IFU

Meta Biomed Co., Ltd-Head Office  
270 Osongsaengmyeong1-ro,  
Osong-eup, Heungdeok-gu,  
Cheongju-si, Chungbuk, Korea

Phone: +82 43 218 1981  
info@meta-biomed.com  
www.meta-biomed.com

Meta Biomed Inc.  
3015 Advance Ln, Colmar,  
Pa 18915, U.S.A

Phone: +1 267 282 5893  
america@metabiomed-inc.com  
www.buymetabiomedonline.com

Meta Biomed Europe GmbH  
Haus der Wirtschaft, Raum 3.11  
Wiesenstraße 35 45473  
Mülheim an der Ruhr, Germany

Phone: +49 208 309 91910  
europe@metabiomed-inc.com  
www.meta-europe.com

**META BIOMED CO.,LTD.**

☎ 043-218-1981  
🌐 www.meta-biomed.com  
📘 www.facebook.com/Metabiomed.Endodontics