

Space-Fill



360° Rotatable Autoclavable Silver Needle

- Innovative coupling technology prevents leakage and overflow
- 999 pure silver maintain thermoplasticity for smooth flow
- Pre-Bent, extra-long design for convenience and visibility



Electric Motorized GP Flow

- Instant start/stop, retraction prevents leakage
- Smooth, bubble-free extrusion for dense & precise filling



Five Memory Settings (T1-T5)

- Quick access to personalized configurations



Precision Temperature Control

- 3 adjustable filling speeds: low, medium, high
- Temperature range: 100°C - 200°C



LED Display Screen


- Real-Time monitoring of gutta-percha remaining





Leakage-Proof Outlet Design

Innovative extrusion outlet structure effectively prevents leakage and simplifies the cleaning process

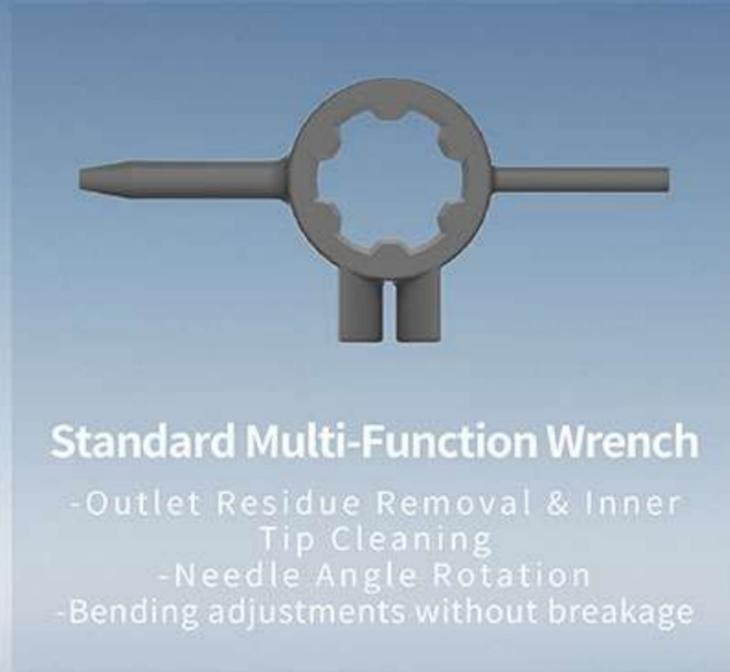

Standard Multi-Function Wrench

- Outlet Residue Removal & Inner Tip Cleaning
- Needle Angle Rotation
- Bending adjustments without breakage

Innovative Two Filling Modes for Diverse Clinical Needs

- Standard Filling Mode
- Quantitative Filling Mode with five selectable settings

Standard Insulating Sleeve for Enhanced Safety

Comprehensive Fill Tip Selection





Magnetic Fluency Consistent Power

Milliohm-level, low resistance for rapid heating and sustained power
Magnetic design enables swift, hassle-free exchange




Dual Battery Configuration for All-Day Productivity


Space-Pack: Powers 800 root canal treatments per battery
Space-Fill: Powers 70 root canal treatments per battery

Lightweight Pen-Style Design

- Slim, ergonomic body for comfortable, effortless grip
- Facilitates easy, precise handling during clinical procedures

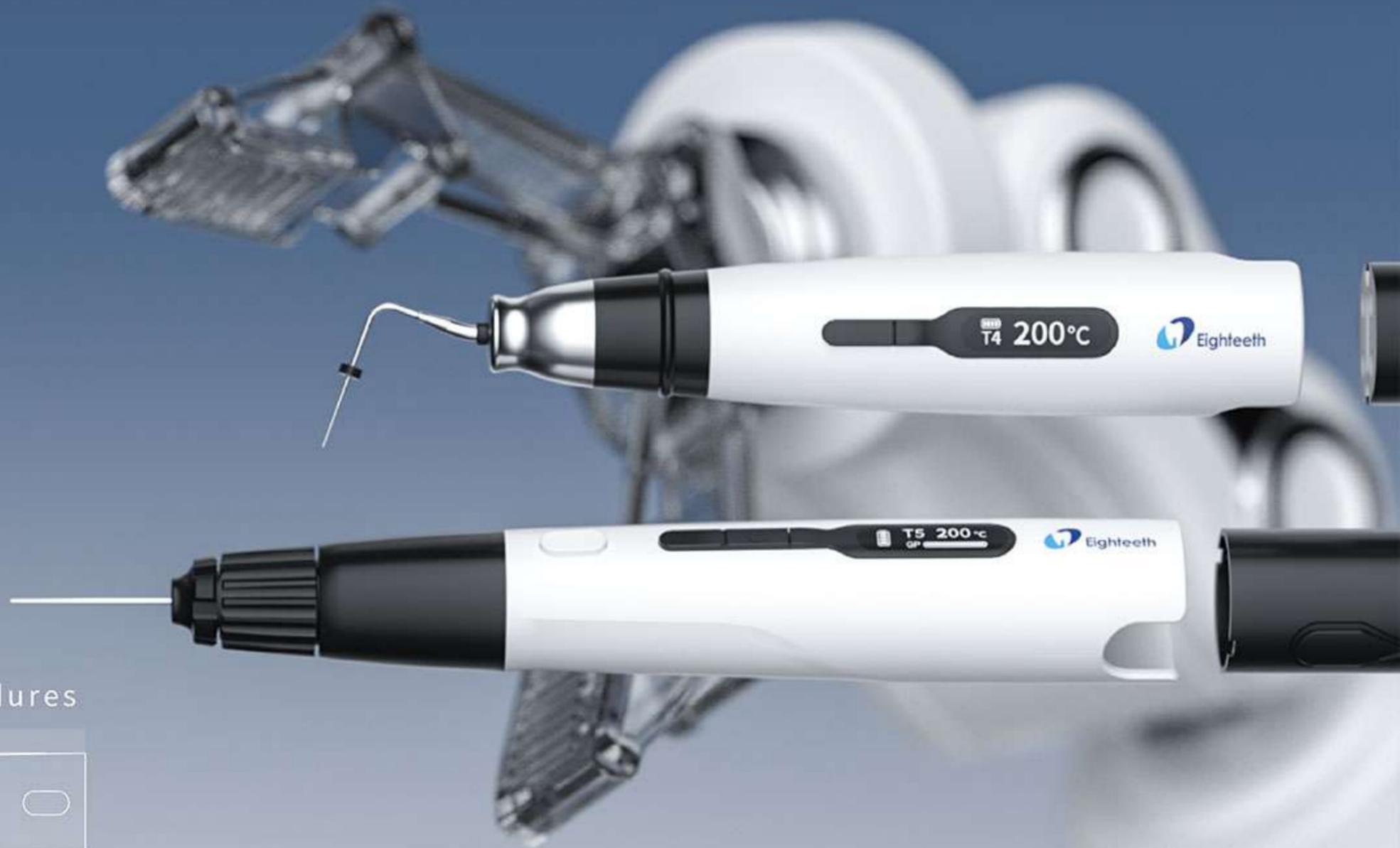

Comprehensive Working Tips for Precision Dentistry

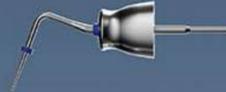Heat Tip			
Standard <i>*(one for each)</i>	Tip 40φ0.4mm	Tip 45φ0.45mm	Tip 50φ0.5mm
			
	Tip 35φ0.35mm	Tip 55φ0.55mm	Tip 60φ0.6mm
Optional			

Fill Tip			
Standard <i>*(one for each)</i>	23Ga	25Ga	
			
Optional			27Ga
			
	23Ga-L	25Ga-L	27Ga-L
			

Specification

Space-Pack	
Heating Time	0.07 seconds for 200 °C (By Tip 35φ0.35mm)
Cooling Time	2 seconds for cooling
Charging Time	45 minutes to complete charge
Operating Time	Up to 800 treatments per full charge
Weight (Include Tip)	84g
Length (Include Tip)	180mm

Space-Fill	
Heating Time	10 seconds for 180°C
Cooling Time	2 seconds for cooling
Charging Time	45 minutes to complete charge
Operating Time	Up to 70 treatments per full charge
Weight (Include Tip)	129g
Length (Include Tip)	206mm

